

ONE-PIECE IMPLANT

FIXO LINE

CATALOGUE 2021





Design & Production	5
Surface Treatment AMS	6
Package	8
General features	10
FIXO MINI diameters and lengths	12
FIXO diameters and lengths	13
FIXO MINI prosthetic components	16
FIXO prosthetic components	17
Surgical instruments	18
FIXO MINI drilling sequence	23
FIXO drilling sequence	24
Guided surgery	26



OXY IMPLANT

DESIGN & PRODUCTION



OXY Implant is a dental implant system totally designed, developed and manufactured in Italy.

The specific company know-how, acquired during thirty years completely spent in researches for dental implantology innovative solutions, allows to keep the **OXY Implant** line devices constantly at the quality level required by the most advanced scientific literature. All this offers the Surgeon a variety of choices that allow to easily deal with even the most complicated clinical situations.

The pursuit of the highest quality, obtained through the careful selection of suppliers and raw materials, the collaboration with research centres and the constant dialogue with the best Italian and foreign Implantologists, results in the great long-term reliability of the **OXY Implant Dental System** with the advantage of the Doctors and Patients ease.

All the implants of the **OXY Implant** line are made of Titanium Grade 4 of European origin, with characteristics of high mechanical resistance resulting from a special cold wire-drawing process.

The prosthetic components are made of Titanium Grade 5 and the surgical instruments of a special stainless steel. The machines used to produce **Oxy Implant** devices are equipped with the best numerical control technology, which allows working with tolerances of a few microns, thus ensuring the excellent overall quality of the implant-abutment system.

The production, testing and packaging phases of all the **OXY Implant** line devices are developed entirely in the company, with the consequent possibility of a direct and constant control of the whole process by a highly skilled Team of Technicians.

The primary packaging of the implants is performed in a cleanroom to avoid any contamination, in compliance with the most stringent sectoral norms.

The renovated head office allows to welcome in an appropriate and modern way all the Customers who wish to observe how the process of realisation of implants, prosthetic components and surgical instruments develops. A large training room allows to organise for Dentists and Dental Technicians many training and refresher events which also represent a useful moment of dialogue between Manufacturer and Users. Their suggestions, derived from daily practice, allow indeed to continuously improve and innovate the **OXY Implant** system.


The medical devices of the **OXY Implant** system comply with the European Directive 93/42/EEC. The manufacturer Biomec S.r.l. has obtained the authorisation to sell from a European notified body and has been equipped with a quality system ISO 9001 and ISO 13485 since 1998.

SURFACE TREATMENT



AMS (Advanced Micro Surface) is the surface selected for the entire range of implants by **OXY Implant**. This surface goes far beyond the results of earlier machined surfaces or surfaces treated with acid and chemical passivation. The **AMS** treatment was developed to accelerate the biological response of cell adhesion, thus enhancing the osseointegration process and the final clinical outcome. **AMS** is obtained with processes of chemical etching, decontamination and a cold Argon plasma treatment. Together they create a new surface with a high degree of cleanliness characterised by a homogeneous micro-porosity of the order of a few microns, smaller than the cell size, so as to considerably increase the adhesion of the osteoblastic cells and to favour the process of osteogenesis.



All these factors are decisive in achieving the best bone-to-implant anchorage and a higher torque required for implant removal. The surface treatment of endosseous components is carried out by a specialised and certified company.

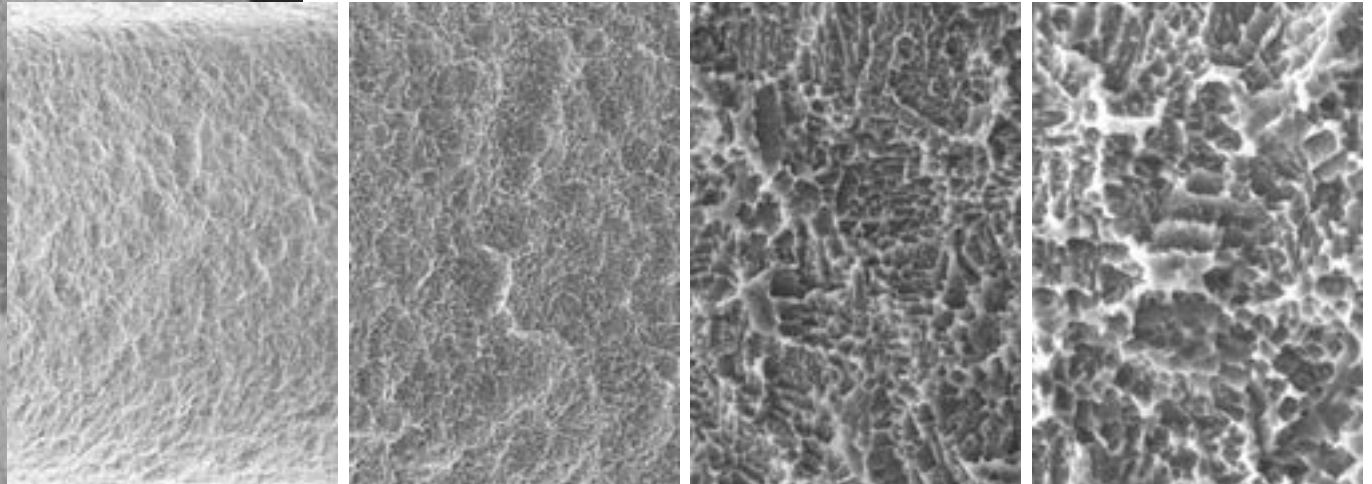

The 5000 X and 7500 X images show in detail the roughness obtained by means of the treatment: the interpeak distance is just a few microns, certainly smaller than the cell size, and accords with recent data about the effects of roughness size on the differentiation and behaviour of osteoblastic cells.

The cleanliness of the screw surface has been confirmed by an XPS analysis on the surface chemical composition. The analysed depth is about 5nm and provides a direct indication of the chemical composition of the material layers directly in contact with the bone.

According to the literature, the maximum concentration of Ti observable by means of XPS on the implant surface after treatment is 14-19%. A percentage of Titanium higher than 10% can be considered satisfactory. The detected amount of Titanium on **OXY** implants is 18%, a value close to the maximum obtainable percentage.

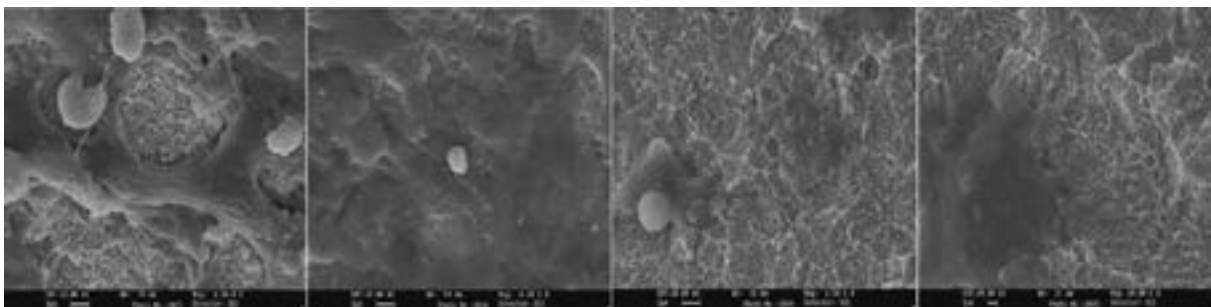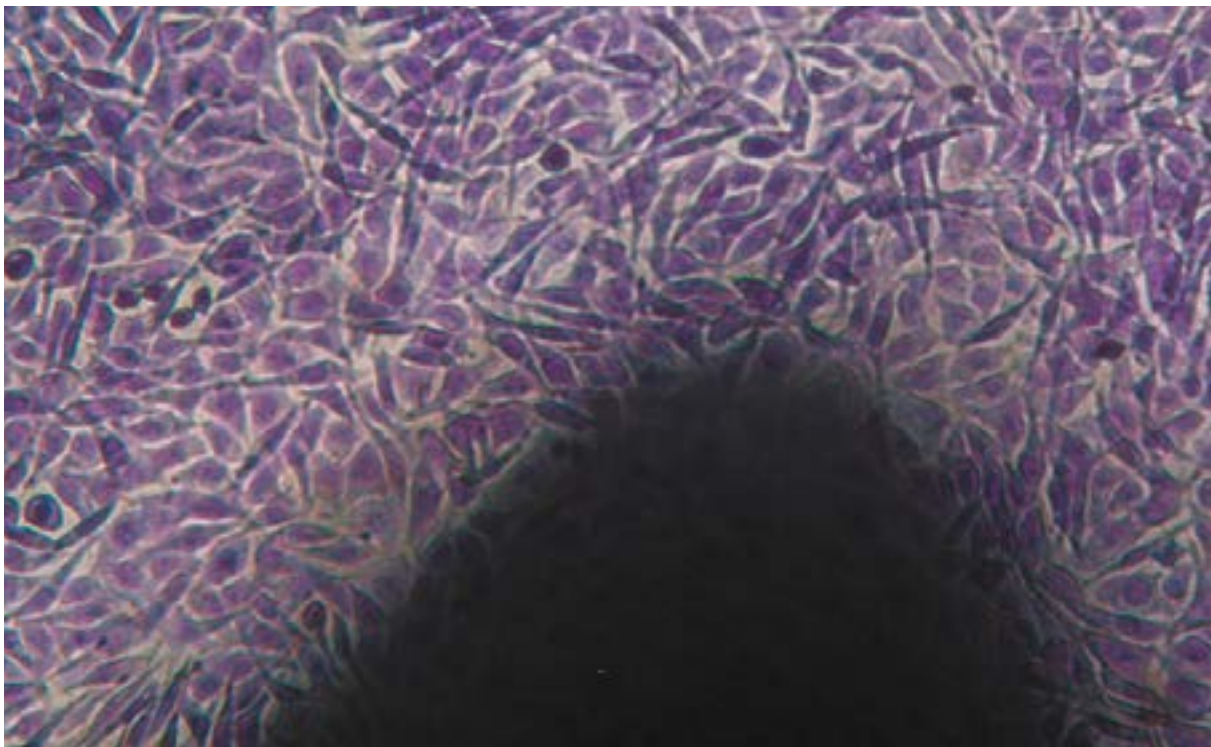

Cytotoxicity tests have been carried out in order to confirm the absence of toxicity induced by the processed screws. The possible presence of dead cells, multinucleated giant cells and general abnormalities in cell morphology was evaluated by comparing the data obtained at the cell layer in contact with the negative control (gold cylinder the same size as the implants) and with the positive control (gutta-percha cylinder).

After being observed under a microscope, the cells were fixed with fixative solution, coloured and photographed. The image relates to the results of tests carried out on experimental samples of the **OXY Implant** line.

All the observations reveal a situation in line with that of the negative control, that is they confirm the absence of toxic effects.

The surface treatment of the **OXY** implants guarantees:

- osteoblast adhesion to the implant surface
- no adverse effect on osteoblast adhesion and growth nor on the surface colonisation.

In conclusion:

- the surface treatment process implemented causes a morphological alteration of the processed screws
- the roughness obtained conforms to what is deemed appropriate to enhance the healing process of bone and osseointegration, at the current state of knowledge
- the implemented process and the following decontamination step allows the complete removal of residues arising from processing and the production of surfaces free of contaminants and foreign deposits
- implants processed with this treatment do not show cytotoxicity effects
- validated and checked periodically, 100% visual inspection and electron microscope analysis on some samples of each treatment batch.

PACKAGE

The implants outer package of **OXY Implant** is made of cardboard and has a convenient tear-off opening system.

The design of the package contains the name of the implant line, the angulation of the Fixo MUA head and the colour coding.

The external stickers of the package indicate all the specific data of the implant: description, reference number, production batch, expiry date and manufacturer.

Inside the package there is the sterile tube containing the implant assembled to the mounter, the instructions for use and 3 stickers indicating all the specific data of the implant, which should be applied respectively on the medical record, on the patient's implant passport and on any additional documentation.



The prosthetic components and the optional or spare surgical instruments are packaged in thermosealed polyethylene bags, with an adhesive label showing all the data of the device: description, reference and batch numbers, colour coding, manufacturer and specific standard symbols for each item.

The prosthetic components and surgical instruments are supplied NON-STERILE and must be sterilised before use.



PRIMARY PACKAGING

All the implants made by **OXY Implant** are contained in a plastic material tube sealed with a screw cap with safety ring, tested to guarantee a 5-year sterility. The sterilisation is performed with BETA rays irradiation validated process, guaranteed and subjected to periodic controls.

MOUNT FREE - OXY CLAMP SYSTEM:

- the mounter preassembled on the implant is housed in a plastic support with titanium core
- allows to see the implant before opening the tube and prevents the movements of the implant facilitating the extraction

Extraction of the implant in 4 simple steps:





Fixo line consists of **single-piece** implants, ideal for the Surgeon who wants to perform a fast immediate load rehabilitation thanks to the **integrated** multi-unit abutment.

Three angulations are available - **0°**, **17°** and **30°** - to allow the surgery and prosthetic team to solve the disparallelisms resulting from the straight and angled implants placement, fully meeting the aesthetic and functional requirements.

The **small** diameter of **Fixo line** implants, plus the perfect **coaxiality** of the fixture-mua monobloc, make them ideal for use in **guided surgery**. Thanks to the characteristics above, there is **no interference** during their insertion through the sleeves of the surgical template and during template removal.

The positioning of the implant with the screw-hole direction of the multi-unit abutment perfectly **corresponding** to the design phase planning is ensured by a mounting system conceived and tested in collaboration with **3Diemme** software-house, in which the stated goal is perfectly centred in the three dimensions by matching precise and evident spatial references obtained through special pre-assembled mounter.

However, the tools available in the dedicated surgical kit allow Fixo to be easily used also in **conventional surgery**: a specific **guide** provides support in creating the implant site with the best angulation, while the special **pins** allow to evaluate previously the most appropriate implant angulation to be used in each single case.

The particular anatomical shape of the neck (pink coloured to favour the aesthetic harmonisation with soft tissues), in addition to the already mentioned small diameter, favour an optimal and fast soft tissues **healing**.

It is evident that, in accordance with a widely consolidated clinical literature, the **absence** of implant-abutment connection cancels the possibility of **micromovements** and consequently helps avoiding bacterial infiltrations.


PRE-ASSEMBLED MOUNTER

The special pre-assembled mounters facilitate the implant placement, even with an angled abutment, allowing to perfectly match the screwing to the axis of the surgical site.

The Ø 5 mm mounter allows the use in **guided surgery** of both **Fixo** and **Fixo Mini** and, thanks to the appropriate template sleeve, also the positioning corresponding to the one planned with the software.

Mounter colour coding:


■ 0° ■ 17° ■ 30°

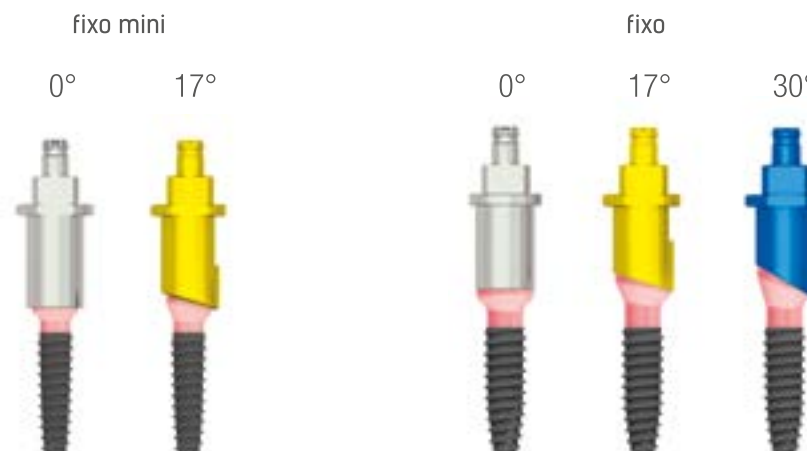

FIXO mini

morphology

- **MORPHOLOGY**
One-piece implant
Chamfer diameter 4.0 mm
Chamfer height: short
- **THREAD**
Single start thread
- **COLOUR CODING**
0° angulation: GREY
17° angulation: YELLOW
- **MATERIAL**
Grade 4 Titanium commercially pure:
medical (ASTM F67)
Cold worked:
resistance more than 40%
Bar tolerance h6:
more precision in components
manufacturing

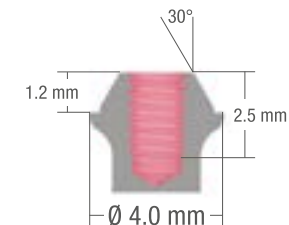


The chamfer diameter of the integrated abutment (0° and 17°) is only **4.0 mm**: Fixo Mini is therefore ideal for use in areas of high **aesthetic** value. Its small diameters allow for use only in low stress incisive areas.

The peculiar shape does not however compromise its mechanical strength and allows the use of a **M1.8** prosthetic screw with tightening to **30 Ncm**.

FIXO

morphology

- **MORPHOLOGY**
One-piece implant
Chamfer diameter 4.8 mm
Chamfer heights: short and long
- **THREAD**
Single start thread
- **COLOUR CODING**
0° angulation: GREY
17° angulation: YELLOW
30° angulation: BLUE
- **MATERIAL**
Grade 4 Titanium commercially pure:
medical (ASTM F67)
Cold worked:
resistance more than 40%
Bar tolerance h6:
more precision in components
manufacturing

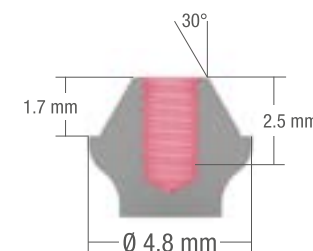

The **prosthetic screw** has an **M1.8 thread** (compared to the standard M1.4 of the two-stage solutions).

This allows the prosthetic components to be tightened with a **torque** of **30 Ncm**.

FIXO mini

diameters and lengths

Ø 3.0 mm

mini



Length	Ref.
10 mm	F4M3000100
11.5 mm	F4M3000115
13 mm	F4M3000130




Length	Ref.
10 mm	F4M3017100
11.5 mm	F4M3017115
13 mm	F4M3017130

Ø 3.5 mm

mini



Length	Ref.
8.5 mm	F4M3500085
10 mm	F4M3500100
11.5 mm	F4M3500115
13 mm	F4M3500130
15 mm	F4M3500150

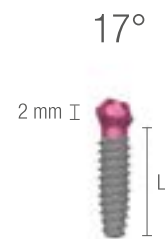

Length	Ref.
8.5 mm	F4M3517085
10 mm	F4M3517100
11.5 mm	F4M3517115
13 mm	F4M3517130
15 mm	F4M3517150


The nominal length L is measured from the apex up to 0,5 mm above the treated section.

FIXO

diameters and lengths

Ø 3.5 mm


short

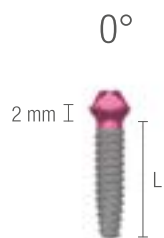


Length	Ref.
8.5 mm	F4G3500085S
10 mm	F4G3500100S
11.5 mm	F4G3500115S
13 mm	F4G3500130S
15 mm	F4G3500150S




Length	Ref.
8.5 mm	F4G3517085S
10 mm	F4G3517100S
11.5 mm	F4G3517115S
13 mm	F4G3517130S
15 mm	F4G3517150S

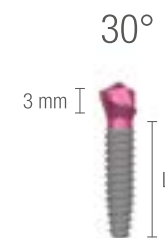


Length	Ref.
10 mm	F4G3530100S
11.5 mm	F4G3530115S
13 mm	F4G3530130S
15 mm	F4G3530150S

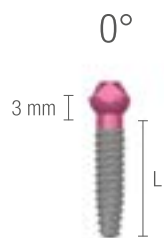

Length	Ref.
10 mm	F4G3500100L
11.5 mm	F4G3500115L
13 mm	F4G3500130L
15 mm	F4G3500150L

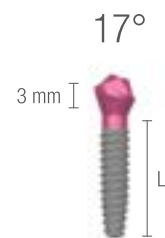

Length	Ref.
10 mm	F4G3517100L
11.5 mm	F4G3517115L
13 mm	F4G3517130L
15 mm	F4G3517150L

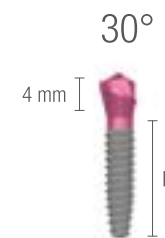

Length	Ref.
10 mm	F4G3530100L
11.5 mm	F4G3530115L
13 mm	F4G3530130L
15 mm	F4G3530150L

The nominal length L is measured from the apex up to 0,5 mm above the treated section.

FIXO


diameters and lengths

Ø 4.0 mm


short




Length	Ref.
7 mm	F4G4000070S
8.5 mm	F4G4000085S
10 mm	F4G4000100S
11.5 mm	F4G4000115S
13 mm	F4G4000130S
15 mm	F4G4000150S
17 mm	F4G4000170S

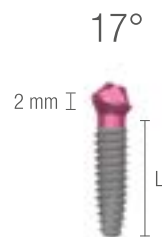


Length	Ref.
7 mm	F4G4017070S
8.5 mm	F4G4017085S
10 mm	F4G4017100S
11.5 mm	F4G4017115S
13 mm	F4G4017130S
15 mm	F4G4017150S
17 mm	F4G4017170S

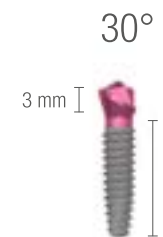


Length	Ref.
10 mm	F4G4030100S
11.5 mm	F4G4030115S
13 mm	F4G4030130S
15 mm	F4G4030150S
17 mm	F4G4030170S

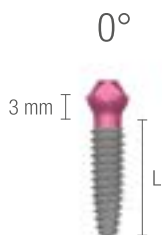

Length	Ref.
10 mm	F4G4000100L
11.5 mm	F4G4000115L
13 mm	F4G4000130L
15 mm	F4G4000150L

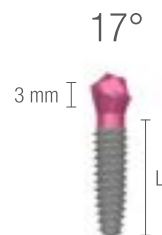

Length	Ref.
10 mm	F4G4017100L
11.5 mm	F4G4017115L
13 mm	F4G4017130L
15 mm	F4G4017150L

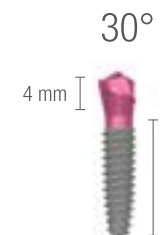

Length	Ref.
10 mm	F4G4030100L
11.5 mm	F4G4030115L
13 mm	F4G4030130L
15 mm	F4G4030150L

The nominal length L is measured from the apex up to 0,5 mm above the treated section.

FIXO

diameters and lengths

Ø 4.5 mm

short

0°



Length	Ref.
7 mm	F4G4500070S
8.5 mm	F4G4500085S
10 mm	F4G4500100S
11.5 mm	F4G4500115S
13 mm	F4G4500130S
15 mm	F4G4500150S

17°



Length	Ref.
7 mm	F4G4517070S
8.5 mm	F4G4517085S
10 mm	F4G4517100S
11.5 mm	F4G4517115S
13 mm	F4G4517130S
15 mm	F4G4517150S

30°



Length	Ref.
10 mm	F4G4530100S
11.5 mm	F4G4530115S
13 mm	F4G4530130S
15 mm	F4G4530150S

Ø 4.5 mm

long

0°



Length	Ref.
10 mm	F4G4500100L
11.5 mm	F4G4500115L
13 mm	F4G4500130L
15 mm	F4G4500150L

17°



Length	Ref.
10 mm	F4G4517100L
11.5 mm	F4G4517115L
13 mm	F4G4517130L
15 mm	F4G4517150L

30°



Length	Ref.
10 mm	F4G4530100L
11.5 mm	F4G4530115L
13 mm	F4G4530130L
15 mm	F4G4530150L

The nominal length L is measured from the apex up to 0,5 mm above the treated section.

FIXO mini

prosthetic components

Prosthetic colour coding
FIXO:

 YELLOW

Implant libraries
available for:

- EXOCAD®
- 3SHAPE®
- DENTAL WINGS®

impression
coping
transfer

Titanium Grade 5
screw included



F4MTR

antirotational
abutment
transfer

Titanium Grade 5
prosthetic screw included



F4MMTA

castable
abutment

PMMA resin
prosthetic screw included



F4MCA

prosthetic
screw
M1.8

Titanio Grado 5
rivestimento black diamond



F4GPS

temporary
abutment

Titanium Grade 5
prosthetic screw included



F4MTA

scan body

Titanium Grade 5



SCANF4M

healing
cap

Titanium Grade 5



F4MHC

prosthetic
screw
with torx

Titanium Grade 5
black diamond coating



F4GCAR

cover
for
bonding

Titanium Grade 5
prosthetic screw included



F4MCT



F4MCC

digital and
laboratory
analogue

Stainless steel



F4MAD

FIXO

prosthetic components

Prosthetic colour coding
FIXO mini:

PINK

Implant libraries
available for:

- EXOCAD®
- 3SHAPE®
- DENTAL WINGS®

impression
coping
transfer



F4GTR

Titanium Grade 5
screw included

antirotational
abutment
transfer



F4GMTA

Titanium Grade 5
prosthetic screw included

castable
abutment



F4GCA

PMMA resin
prosthetic screw included

prosthetic
screw
M1.8



F4GPS

Titanio Grado 5
black diamond coating

temporary
abutment



F4GTA

Titanium Grade 5
prosthetic screw included

scan body



SCANF4G

Titanium Grade 5

healing
cap



F4GHC

Titanium Grade 5

prosthetic
screw
with torx



F4GCAR

Titanium Grade 5
black diamond coating

cover
for
bonding



F4GCT



F4GCC

Titanium Grade 5
prosthetic screw included

oxy loc
abutment



F4GLOC

Titanium Grade 5
Gold Titanium coating

Locator® compatible
Locator® is a Zest Anchors
registered trademark

digital and
laboratory
analogue



F4GAD

Stainless steel

mechanical
screwdriver
for angled
screwing



SESD



MESD

torx key

Stainless steel
black diamond coating

Mechanical screwdriver SESD or MESD must be ordered separately.

SURGICAL INSTRUMENTS



The drills for **Fixo** and **Fixo mini** implants are designed and manufactured with the aim of ensuring the constant efficiency and durability. These prerogatives are realised firstly thanks to the use of absolute quality materials as the special hardened stainless steel used for their manufacture. Finishing processes of last generation complete the production cycle.

All the drills are **Black Diamond** coated, a perfectly biocompatible material which increases their superficial hardness and decreases the generated friction forces. This process increases therefore considerably the cutting power and consequently reduces the risk of the bone overheating during surgery.

Their morphology is cylindrical with a smaller diameter step to provide a surgical socket that follows accurately the shape of the **Fixo** and **Fixo mini** implants, more conical in the apical area.

The drills are sequential with a 0.3 mm step between the successive diameters (2.9 - 3.2 - 3.5 - 3.8 and 4.1 mm) to create a surgical socket corresponding to the diameter of the implant core or to underprepare or overprepare it according of the bone density characteristics.

In order to facilitate the Surgeon's work, they are graduated with laser markings to indicate the length of the corresponding implant to be inserted. For the same reason they can also be provided with titanium depth stoppers.

Also, the diameter of the drill is well highlighted thanks to a specific laser marking.

The colour of the ring on the shank of the instrument allows to easily locate the housing of the surgical kit in which to place it: it is in fact characterised by the same coding.

It is important to remember that the tip of the drill increases of 0.5 mm the nominal length of the instrument. Taking always this into account, will be possible to avoid the anatomical structures damage especially during use in areas close to the maxillary sinus or to the mandibular channel.



SURGICAL INSTRUMENTS

pilot drill

Ø 1.8 mm

Hardened stainless steel
black diamond coating

max 900 Rpm



SID01L
long

initial graduated drill

Ø 1.6 mm

Hardened stainless steel
black diamond coating

max 1000 Rpm



3DM006131615

initial graduated drill

Ø 2.0 mm

Hardened stainless steel
black diamond coating

max 1000 Rpm



IGSD200

initial graduated drill

Ø 2.3 mm

Hardened stainless steel
black diamond coating

colour coding: white

max 900 Rpm



IGSD230

final graduated drill

Ø 2.9 mm

Hardened stainless steel
black diamond coating

colour coding: yellow

max 900 Rpm



SD2923

final graduated drill

Ø 3.2 mm

Hardened stainless steel
black diamond coating

colour coding: red

max 800 Rpm



SD3226

final graduated drill

Ø 3.5 mm

Hardened stainless steel
black diamond coating

colour coding: green

max 800 Rpm



SD3529

final graduated drill

Ø 3.8 mm

Hardened stainless steel
black diamond coating

colour coding: blue

max 700 Rpm



SD3832

final graduated drill

Ø 4.1 mm

Hardened stainless steel
black diamond coating

colour coding: brown

max 700 Rpm



SD4135

calibrated expandable stop

Titanium Grade 5

Length	Ref.
5.5 mm	EST4055T
7 mm	EST4070T
8.5 mm	EST4085T
10 mm	EST4100T
11.5 mm	EST4115T
13 mm	EST4130T
15 mm	EST4150T



parallel indicator

Titanium Grade 5



PIN2329
ø2.3-2.9 mm

PIN2317
17°

PIN2330
30°

surgical direction guide

Anodised Titanium



ILSGUIDE

countersink drill

Hardened stainless steel
black diamond coating
max 300 Rpm



STDPSILS

tissue punch

Hardened stainless steel
black diamond coating
max 100 Rpm



MPU34

driver converter square 4x4

grey disc
Hardened stainless steel



XCMRC8
extrashort

CMRC8
short

LMRC8
medium

drill extension

Hardened stainless steel



DEXT

mechanical adapter fixo

Hardened stainless steel
Tightening torque 50 Ncm



MCF4G

torque ratchet

Stainless steel
driver D8
torque 10÷70 Ncm



TR8

SURGICAL INSTRUMENTS

The surgical kits are designed with the purpose to enhance the compactness and ease of use.

The arrangement of the instruments and the colour coding determines a logical and intuitive instruments use process during the successive steps of the intervention.

To facilitate their maintaining, the boxes are made of autoclavable and autosanitising plastic.



ref. TRF4G

additional surgical kit:

ref. TREIPSK

ref. TREIPSCLS



ref. TRPSKFIXO


surgical instrument kit

for \varnothing 3.0 mm FIXO MINI implants

D1-D2 BONE

pilot
 \varnothing 1.8 mm
900 Rpm

graduated
 \varnothing 2.0 mm
1.000 Rpm

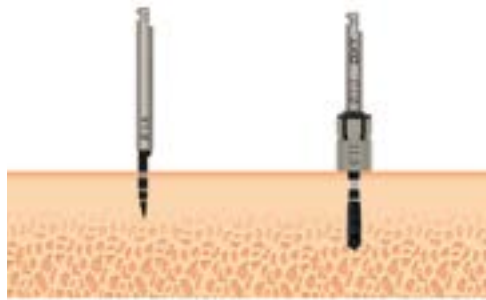

D3-D4 BONE

pilot
 \varnothing 1.8 mm
900 Rpm

graduated
 \varnothing 1.6 mm
1.000 Rpm



Important: the \varnothing 2.0 mm drill is included in the kit, and the \varnothing 1.6 mm drill (ref. 3DM006131615) must be purchased separately.

for \varnothing 3.5 mm FIXO MINI implants

pilot
 \varnothing 1.8 mm
900 Rpm

initial graduated
 \varnothing 2.3 mm
900 Rpm

parallel
indicator
 \varnothing 2.3 mm

final graduated
 \varnothing 2.9 mm + stop L 7
900 Rpm


final graduated
 \varnothing 2.9 mm
900 Rpm




for D3-D4 bone

for D1-D2 bone

positioning the ILS GUIDE

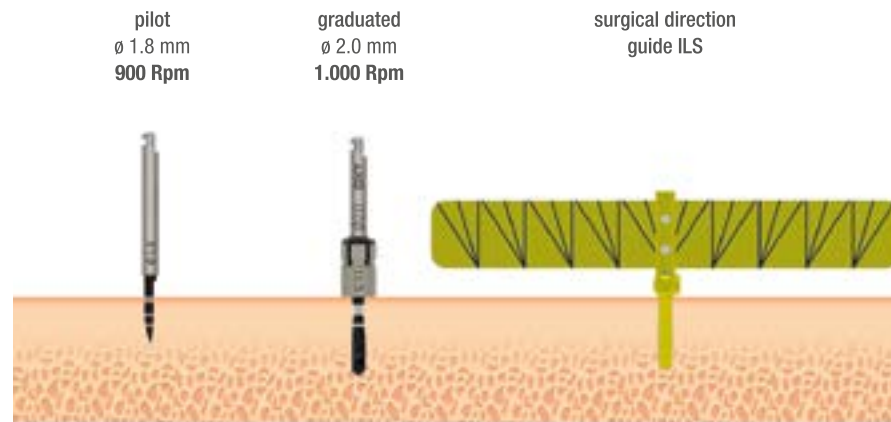


for ø 3.5 mm FIXO implants

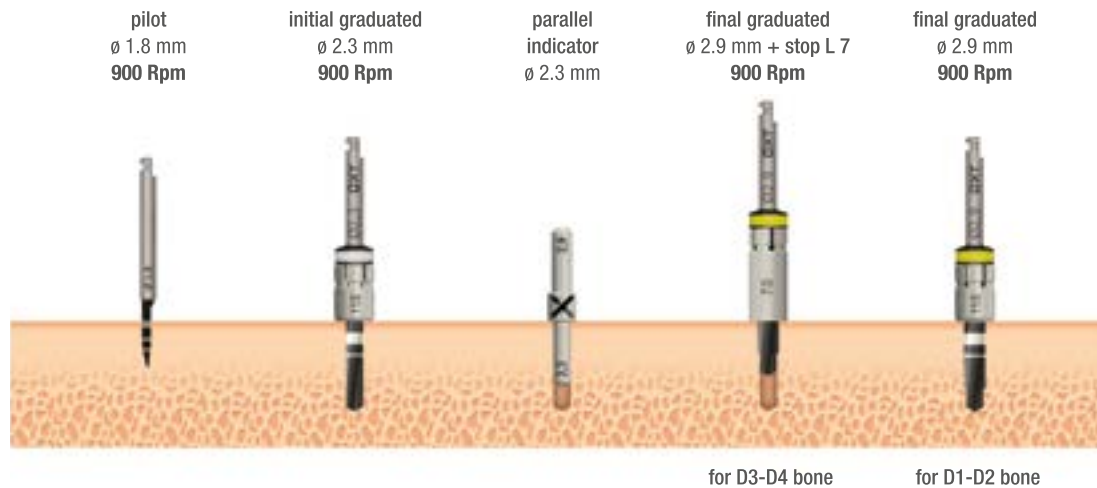


FIXO

DRILLING SEQUENCE


for \varnothing 4.0 mm FIXO implants

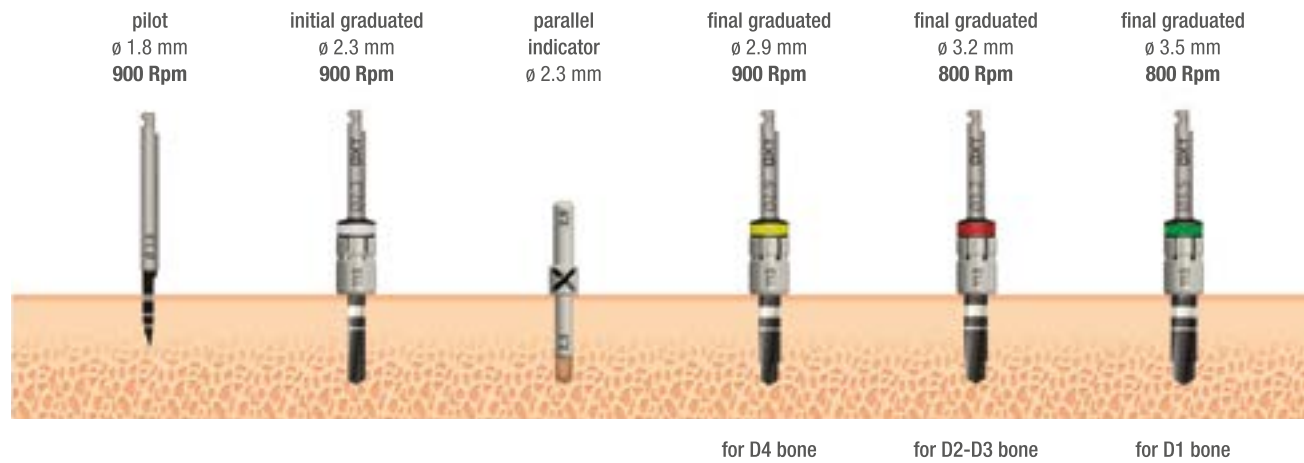

for \varnothing 4.5 mm FIXO implants

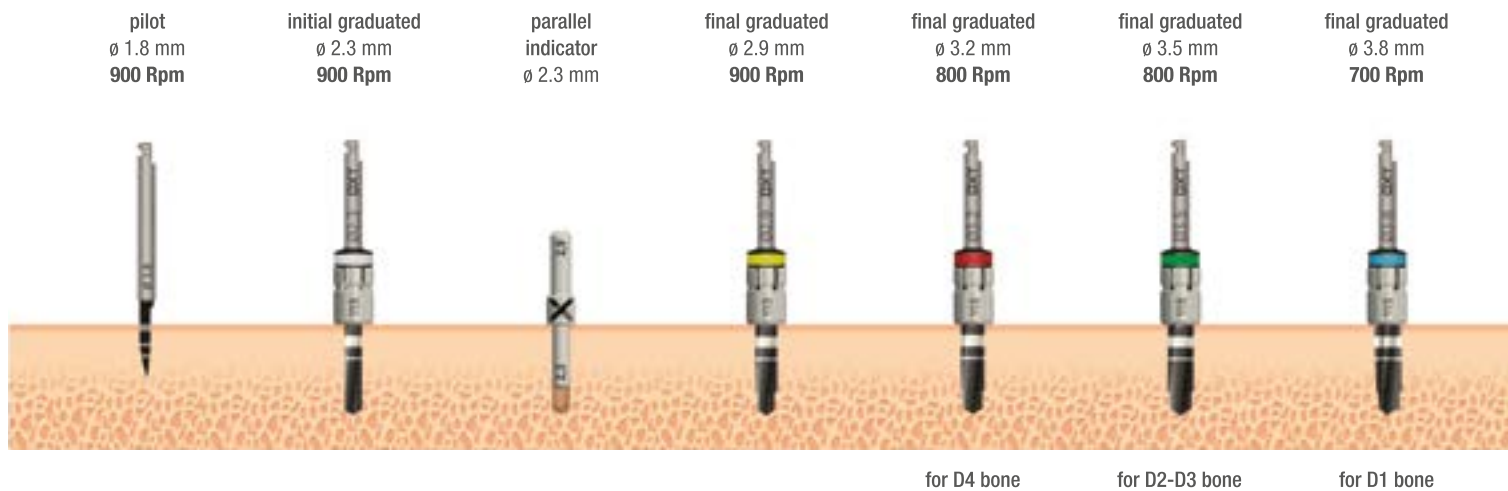

GUIDED SURGERY

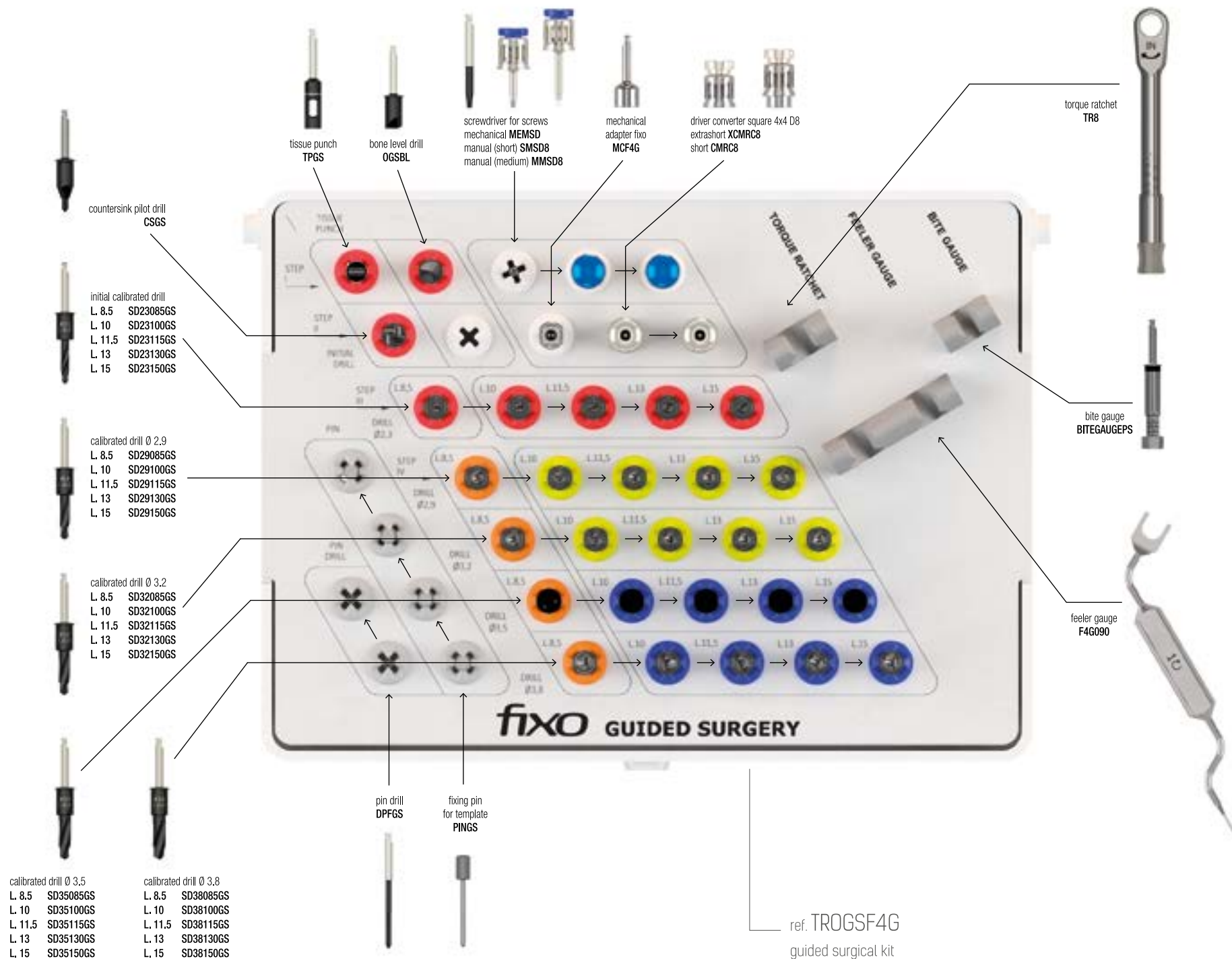

pin drill

Hardened stainless steel
black diamond coating
max 1000 Rpm



DPFGS

tissue punch

Hardened stainless steel
black diamond coating
max 100 Rpm



TPGS

bone level drill

Hardened stainless steel
black diamond coating
max 300 Rpm



OGSBL

countersink pilot drill

Hardened stainless steel
black diamond coating
max 600 Rpm



CSGS

initial calibrated drill

Ø 1.6 mm

Hardened stainless steel
black diamond coating
max 1000 Rpm



Length	Ref.
7 mm	3DM006131670
11.5 mm	3DM0061316115
15 mm	3DM006131615

initial calibrated drill

Ø 2.3 mm

Hardened stainless steel
black diamond coating
max 900 Rpm



Length	Ref.
6 mm	SD23060GS
8.5 mm	SD23085GS
10 mm	SD23100GS
11.5 mm	SD23115GS
13 mm	SD23130GS
15 mm	SD23150GS

calibrated drill

Ø 2.9 mm

Hardened stainless steel
black diamond coating
max 900 Rpm



Length	Ref.
8.5 mm	SD29085GS
10 mm	SD29100GS
11.5 mm	SD29115GS
13 mm	SD29130GS
15 mm	SD29150GS

calibrated drill

Ø 3.2 mm

Hardened stainless steel
black diamond coating
max 800 Rpm



Length	Ref.
8.5 mm	SD32085GS
10 mm	SD32100GS
11.5 mm	SD32115GS
13 mm	SD32130GS
15 mm	SD32150GS

calibrated drill

Ø 3.5 mm

Hardened stainless steel
black diamond coating
max 800 Rpm



Length	Ref.
8.5 mm	SD35085GS
10 mm	SD35100GS
11.5 mm	SD35115GS
13 mm	SD35130GS
15 mm	SD35150GS

calibrated drill

Ø 3.8 mm

Hardened stainless steel
black diamond coating
max 700 Rpm



Length	Ref.
8.5 mm	SD38085GS
10 mm	SD38100GS
11.5 mm	SD38115GS
13 mm	SD38130GS
15 mm	SD38150GS

fixing pin for template

Stainless steel



PINGS

mechanical adapter fixo

Hardened stainless steel
maximum
tightening torque 50 Ncm



MCF4G

driver converter square 4x4

grey disc
Hardened stainless steel



XCMRC8
extrashort



CMRC8
short

mechanical screwdriver for screws

Hardened stainless steel
with Black Diamond coating
Hexagonal key 1.25 mm



MEMSD

manual screwdriver for screws

blue disc
Hardened stainless steel
hexagonal key 1.25 mm



SMSD8
extrashort



CMRC8
short

sleeve for template

Stainless steel



F4GS

fixing pin sleeve

Stainless steel



CPINGS

bite gauge

Stainless steel



BITEGAUGEPS

torque ratchet

Stainless steel

driver D8

torque 10÷70 Ncm



TR8

feeler gauge

Stainless steel



F4G090

NOTES

